


Kyri Cart Evo

MOBILE INSTRUMENT HOLDER UNIT

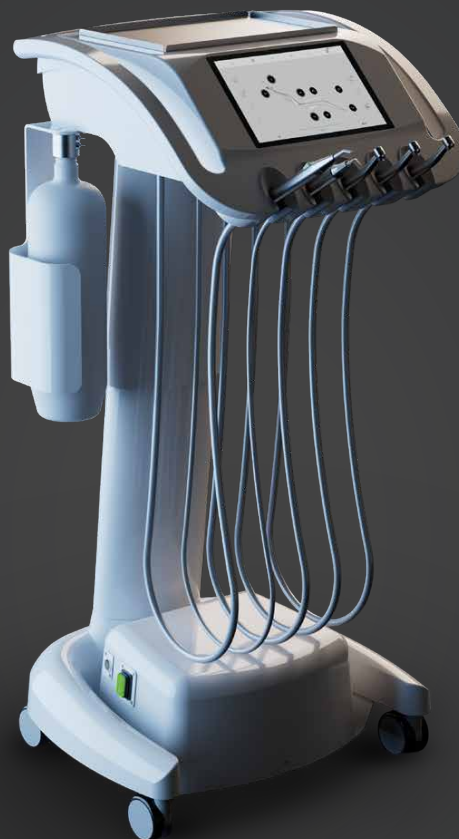

Maximum flexibility, total connectivity.

KYRI EVO AND KYRI CART EVO.

Full integration with the Kyri Evo chair.

Simplify your work and offer
a more comfortable experience for
your patients with a perfectly
connected treatment station.

curved and elegant lines

new design

compact dimensions,
small footprint.


instrument cooling

electronic cooling

with proportional valve control
and real-time display.

KYRI EVO AND KYRI CART EVO. THE INTEGRATED SOLUTION.

It adapts to you. Always.

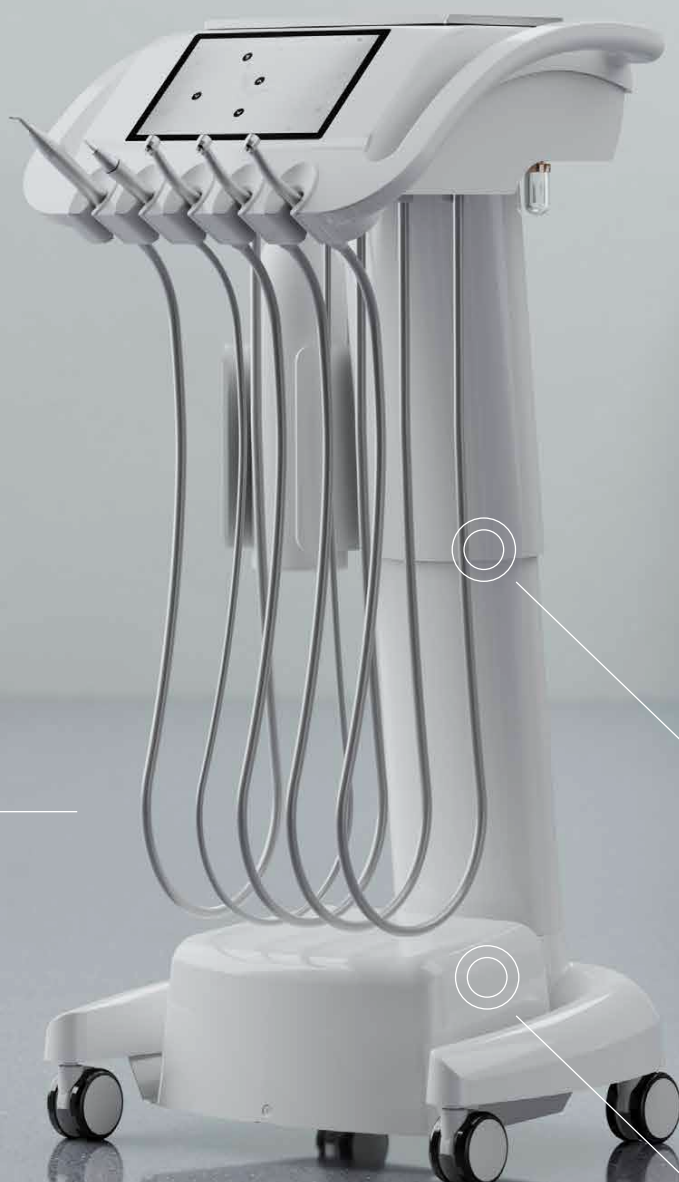

metal protection

new bottle

with metal coupling and protective casing for greater strength, safety and reliability in everyday use.



adaptation to the user profile

variable height from 83 to 103 cm

improved ergonomics,
ease of parking.

connectivity

wi-fi integrated

superior wireless performance for
stable and immediate connectivity.

NEW USER INTERFACE


Intuitive and aesthetic interface.

Every professional can adapt the system to their own operating style: transmission ratios, speed, torque and spray are fully configurable according to the specific requirements of each procedure.

smart interface

new 9" colour touchscreen display


High-resolution display that improves visibility and control

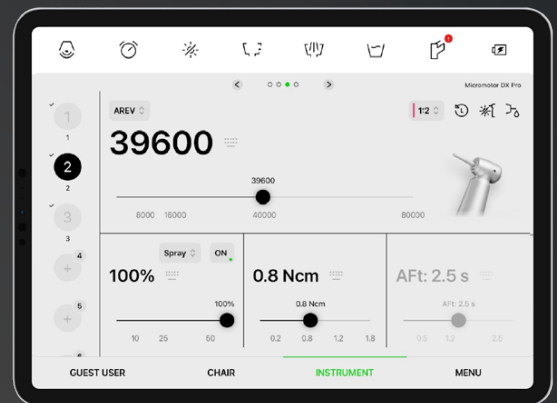

greater versatility

unlimited user profiles

Ability to switch easily from one user to another

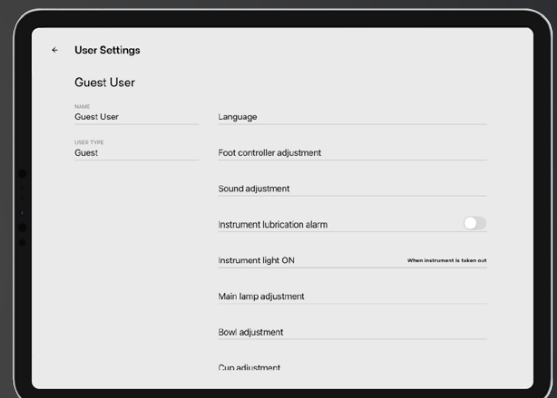


Integration with Kyri Evo

Smart pedal Wireless

Improve your ergonomics by using hands-free operation. Wireless or wired. Single pedal Kyri Cart Evo / Kyri Evo

** The NOK pedal is available.

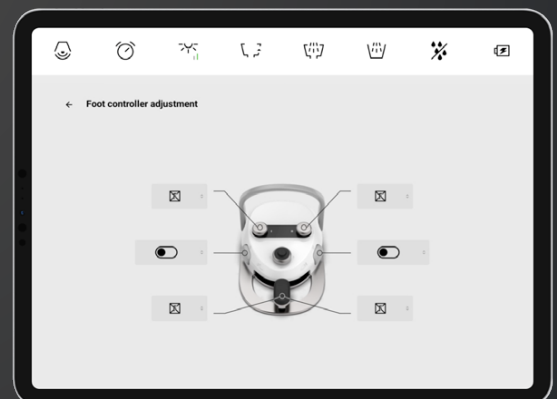

INTERCONNECTION

Maximum flexibility, total connectivity.

Variable connection options:

- Direct connection to the Kyri Evo chair
- Quick connection box
- Ambidextrous connection option



high hygiene standards

function hygiene

Designed to optimise and simplify the cleaning of instrument tubes.

additional slots


three slots for micromotors

Quick access to multiple instruments.
Option to add a surgical micromotor.
Or to combine 2 micromotors and 1 surgical micromotor.

maximum flexibility


operational mobility

Standalone or connected to the central system, it ensures advanced functionality and maximum mobility in any dental environment.

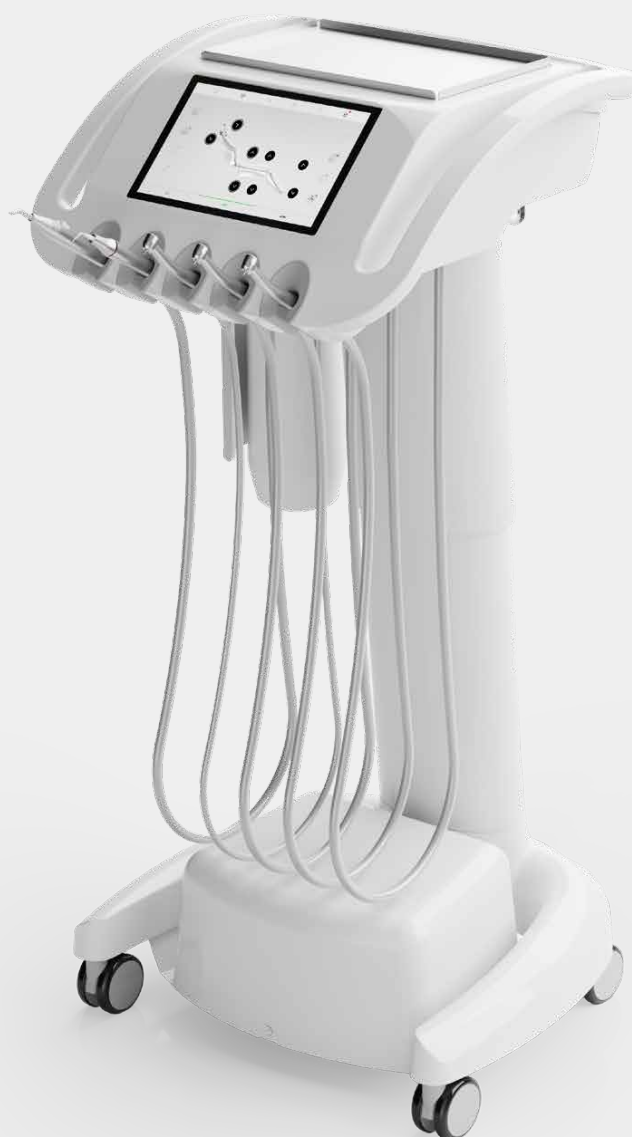


THE VERSATILE SOLUTION

Endless configuration possibilities

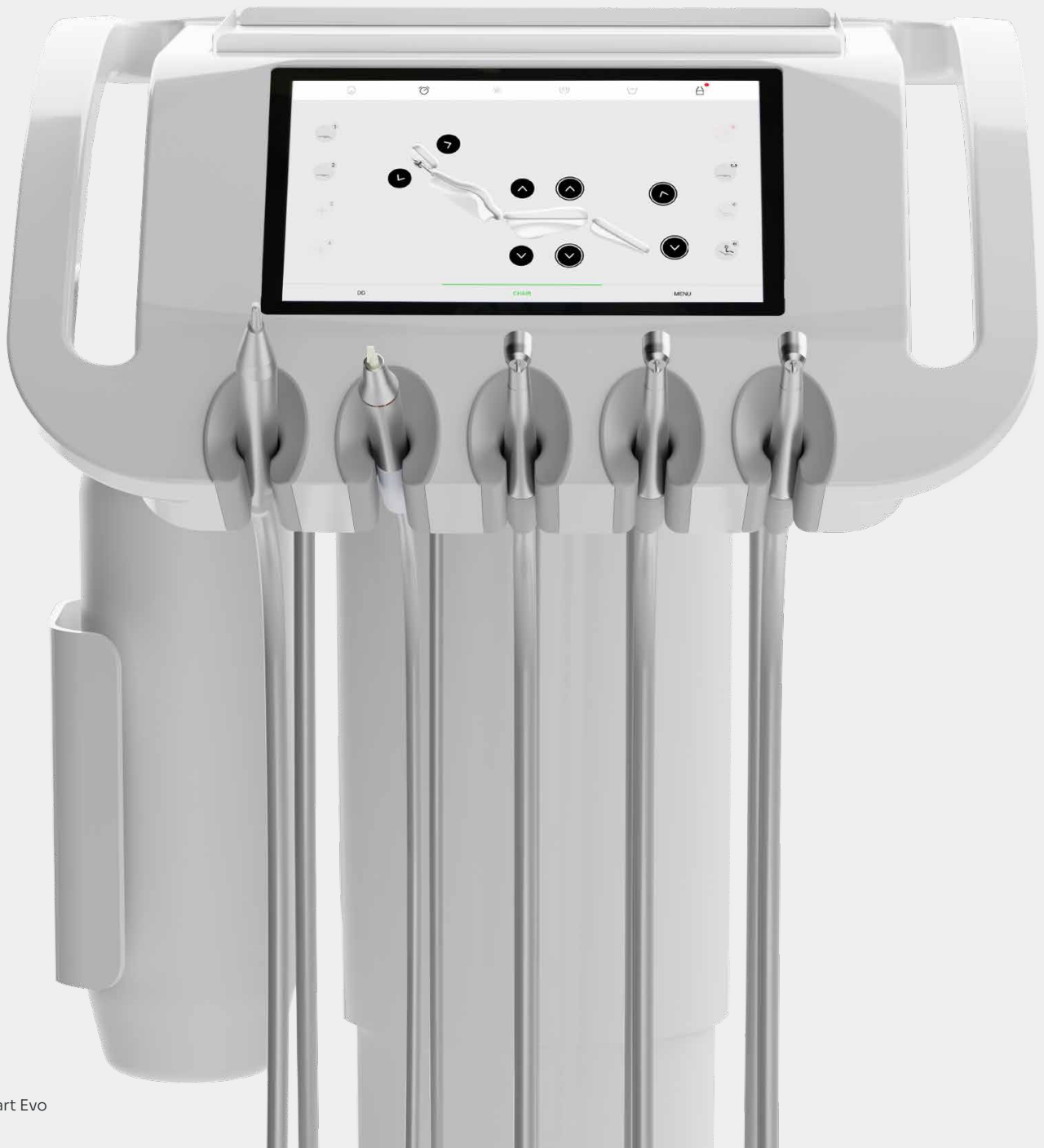


TECHNICAL SPECIFICATIONS

Kyri Cart Evo

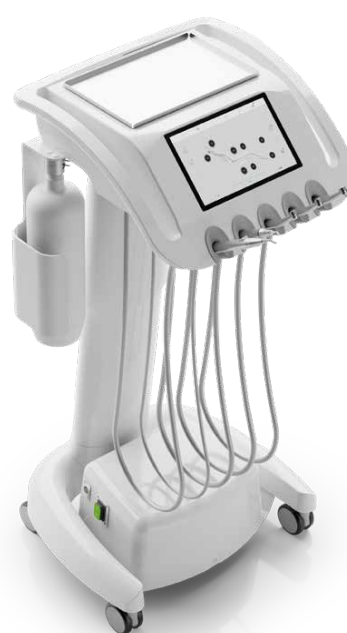

Min./max. height	630 mm – 830 mm
Frequency	50/60 Hz
Maximum input	400 VA
Inlet air pressure	0.45-0.8 Mpa
Inlet water pressure	0.3-0.6 Mpa
Net weight	43 kg (+ max 5kg.)
Width	500 mm
Depth	495 mm
Power supply	220-240 V
Electrical shock protection rating	Protection class I
Cover protection rating	IP21

Kyri Cart Evo

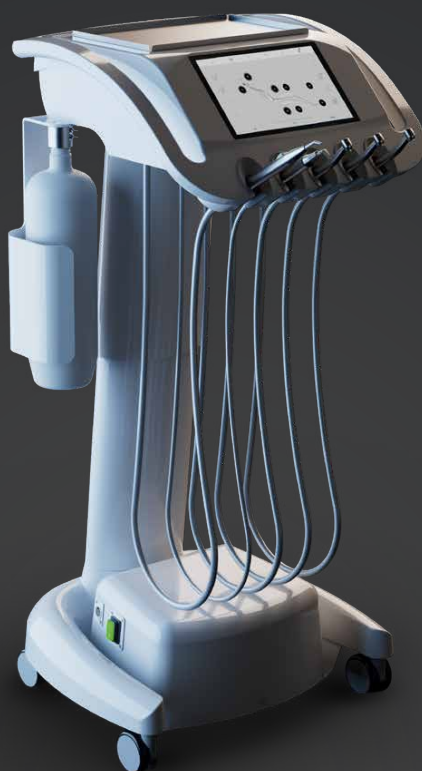

MOBILE INSTRUMENT HOLDER UNIT

Maximum flexibility, total connectivity.



Tecno-Gaz S.p.A.

Strada Cavalli, 4 - 43038 - Sala Baganza - Parma - Italia
Tel. +39 0521 83.80 Fax +39 0521 83.33.91 - www.tecnogaz.com

All rights reserved to Tecno-Gaz S.p.A. Changes to images or content may be made without prior notice.
Tecno-Gaz S.p.A., shall not be considered liable for damage resulting from the lack or inaccuracy of the
information contained in this document. The images used in this document are for illustrative purposes only.

